

O₂ Gas Sensor

(Order Code O2-BTA)



The O₂ Gas Sensor can be used to monitor gaseous oxygen levels in a variety of biology and chemistry experiments. The sensor is intended for measuring gaseous, not aqueous, O₂ concentrations.

Note: Vernier products are designed for educational use. Our products are not designed nor are they recommended for any industrial, medical, or commercial process such as life support, patient diagnosis, control of a manufacturing process, or industrial testing of any kind.

What's Included

- O₂ Gas Sensor
- 250 mL gas sampling bottle (Nalgene bottle with lid)

Compatible Software and Interfaces

See www.vernier.com/manuals/o2-bta for a list of interfaces and software compatible with the O₂ Gas Sensor.

Getting Started

1. Connect the sensor to the interface (LabQuest Mini, LabQuest 2, etc.).
2. Start the appropriate data-collection software (Logger *Pro*, Logger Lite, LabQuest App) if not already running, and choose New from File menu.

The software will identify the sensor and load a default data-collection setup. You are now ready to collect data.

If you are collecting data using a Chromebook™, mobile device such as iPad® or Android™ tablet, or a Vernier wireless sensor or interface, please see the following link for up-to-date connection information:

www.vernier.com/start/o2-bta

Using the Product

Connect the sensor following the steps in the Getting Started section of this user manual. **Note:** For best results, keep the sensor upright as much as possible during use. Reading may drift if used in a horizontal position.

Videos

View videos related to this product at www.vernier.com/o2-bta

Calibration

For many measurements, it will not be necessary to calibrate the O₂ Gas Sensor. The sensor is shipped with a stored calibration.

For more accurate measurements, the sensor can be calibrated at 0 and 20.9% oxygen. Follow the 2-point calibration procedure. For the first point, push and hold the zero button with a pointed object such as a straightened paper clip.

The button is located on the side of the sensor next to CAL. Enter a value of 0 for this reading. Release the button and take a second reading. Enter a value of 20.9% oxygen or a corrected value from the table in the Effect of Humidity section of this user manual. Once finished, the sensor should now read 20.9% (or the value entered from the table). For detailed instructions on calibration see the related links below.

You can also store your custom calibration on the sensor using the following steps:

1. Go to the Storage tab and click the Set Sensor Calibration button (Logger *Pro*) or Save Calibration to Sensor (LabQuest). Once saved, this new calibration will be used automatically, regardless of the interface to which the sensor is connected.
2. Click Done to complete the calibration process.

- For instructions on calibration using Logger *Pro* computer software, see www.vernier.com/til/2341

- For instructions on calibration using LabQuest App, see www.vernier.com/til/3394

- For instructions on calibration using Graphical Analysis with a Chromebook, see www.vernier.com/til/3631

- For instructions on calibration using Graphical Analysis with an iOS or Android device, see www.vernier.com/til/3630

Effect of Humidity

Because the oxygen concentration varies with the amount of water vapor in the atmosphere, you may want to adjust your atmospheric oxygen calibration value to improve accuracy when using the O₂ Gas Sensor. The accepted value of 20.9% for atmospheric oxygen levels is calculated in dry air (0% humidity). If you know the relative humidity of the location at which you are calibrating, you can substitute one of the values below in place of 20.9%.

Relative humidity	0%	25%	50%	75%	100%
Oxygen in % by volume	20.9	20.7	20.5	20.3	20.1

Specifications

Range	0–27% (0–270 ppt)
Accuracy	± 1% volume O ₂ (at 760 mmHg)
Resolution	0.01%
Warm up time	< 5 seconds
Response time	90% of final reading in ~12 seconds
Output signal range	0 to 4.8 VDC; 2.7 to 3.8 VDC @ 21% O ₂
Output impedance	1 KΩ
Input voltage	5 VDC +/-0.25 VDC
Gas sampling mode	Diffusion
Normal operating temperature range	20 to 30°C
Storage temperature range	-20 to +60°C
Pressure range	0.5 to 1.5 atm
Pressure effect	Directly proportional $V_{out} = V_{out}(\text{standard}) \times (P/1013)$
Default calibration values	slope: 6.5625 %/V intercept: 0%
Dimensions	Sensor tube: 45 mm length, 28 mm OD Total length: 120 mm

Care and Maintenance

Important: The O₂ Gas Sensor must be stored upright when not in use. This is necessary to maintain the sensor. Failure to store upright will shorten the life of the sensor.

As your O₂ Gas Sensor ages, the readings will decrease. This is normal, as the chemicals in the electrochemical cell are depleted. It does not mean the sensor is no longer functional; rather, it simply requires that you perform a calibration and store it as described previously. Once the readings in air are below 12% prior to calibration, or 1.8 volts during calibration, the sensor should be replaced. You may choose to have your sensor refurbished and a new element installed, or you may purchase a new sensor.

For tips on how to extend the life of the O₂ Gas Sensor, see www.vernier.com/til/2588

How the Sensor Works

The Vernier O₂ Gas Sensor measures the oxygen concentration in the range of 0

to 27% using an electrochemical cell. The cell contains a lead anode and a gold cathode immersed in an electrolyte. Oxygen molecules entering the cell are electrochemically reduced at the gold cathode. This electrochemical reaction generates a current that is proportional to the oxygen concentration between the electrodes. The sensor output is a conditioned voltage proportional to the reaction current.

Troubleshooting

- Even though the sensor responds rather quickly to changes in O₂ concentration, the gas has to diffuse into the electrochemical cell located at the top of the sensor shaft before any changes in concentration can be detected. Because the diffusion of a gas is a slow process, there can be some delay in readings.
- To collect data in a controlled environment, we recommend that you use the 250 mL Nalgene collection bottle that is included with your sensor. Place the tip of the sensor into the opening of the bottle and push the sensor into the bottle. When the sensor will go no further, you have a seal. **Very important:** Do not place the sensor into any liquid. The sensor is intended only for measuring gaseous, not aqueous, O₂ concentration.
- To collect data in a controlled environment and also use our CO₂ Gas Sensor, we recommend purchasing either the BioChamber 250 (250 mL capacity) or the BioChamber 2000 (2 L capacity). Each BioChamber has a second, grommetted opening so that you can insert two probes at once.

Repair Information

If you have watched the related product video(s), followed the troubleshooting steps, and are still having trouble with your O₂ Gas Sensor, contact Vernier Technical Support at support@vernier.com or call 888-837-6437. Support specialists will work with you to determine if the unit needs to be sent in for repair. At that time, a Return Merchandise Authorization (RMA) number will be issued and instructions will be communicated on how to return the unit for repair.

For information on common repairs for the O₂ Gas Sensor, see www.vernier.com/til/1605

Accessories/Replacements

Item	Order Code
BioChamber 250	BC-250
BioChamber 2000	BC-2000
250 mL Nalgene Bottle w/Lid	CO2-BTL

The BioChambers have two openings with grommets so that you can insert both the CO₂ Gas Sensor and the O₂ Gas Sensor to take measurements of the atmosphere inside the BioChamber simultaneously.


Warranty

Vernier warrants this product to be free from defects in materials and workmanship for a period of five years from the date of shipment to the customer. This warranty does not cover damage to the product caused by abuse or improper use. This warranty covers educational institutions only.

Rev. 6/9/16

Logger *Pro*, Logger *Lite*, Vernier LabQuest 2, LabQuest Mini, and other marks shown are our trademarks or registered trademarks in the United States.

All other marks not owned by us that appear herein are the property of their respective owners, who may or may not be affiliated with, connected to, or sponsored by us.

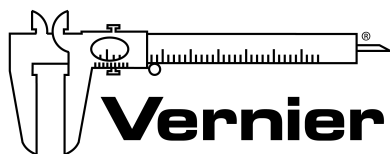

Vernier Software & Technology
13979 SW Millikan Way • Beaverton, OR 97005-2886
Toll Free (888) 837-6437 • (503) 277-2299 • Fax (503) 277-2440
info@vernier.com • www.vernier.com